ENGLISH

Sonicare For Kids (Fig. 1)

- A Hygienic travel cap
- B Brush head
- C Handle with soft grip
- D Power on/off button with battery charge indicator
- E Interchangeable panel
- F Brushing mode button with mode lights
- G Charger with brush head holder

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USE

DANGERS

To reduce the risk of electrocution:

- Do not place or store the product where it can fall or be pulled into a tub or sink.
- Do not place the product in or drop it into water or any other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not use the product while bathing.

WARNINGS

To reduce the risk of burns, electrocution, fire to persons:

- Discontinue use if this product appears damaged in any way (brush head, handle, charger). This product contains no serviceable parts. Refer to chapter 'Guarantee and service' if your Sonicare no longer works properly or needs repair.
- Never use the charger if it has a damaged cord or plug.
- Keep the cord away from heated surfaces.
- Do not use the charger outdoors.
- Do not use the charger if dropped into water.
- This product is designed to clean your teeth, gums and tongue only. Use this product only for its intended use as described in this booklet. Discontinue use of this product and contact a physician/dentist if discomfort or pain is experienced.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not use attachments other than those recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean the brush head and the handle with soap and water after each use. Otherwise the plastic may crack.

MEDICAL WARNINGS

- Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.
- Contact your dental professional if excessive bleeding occurs after using this product or bleeding continues to occur after 1 week of use.
- If you have a pacemaker or other implanted device contact your physician or the device manufacturer prior to use.
- Sonicare has been tested and is compliant with safety standards for electromagnetic devices.
- Consult your physician prior to using the Sonicare if you have medical concerns.

NOT A TOY

CHECK BRUSH HEADS REGULARLY, CRACKS FROM DROPPING CAN LEAD TO CHOKE HAZARD

SAVE THESE INSTRUCTIONS

Getting started

Applying the interchangeable panel

The Sonicare For Kids comes with 3 different interchangeable panels.

- Remove the adhesive backing from the panel.
- Apply the panel to the white section on the front of the handle (Fig. 2).

Note: Make sure you line up the cut-out areas on the panel with the buttons on the front of the handle.

Reapplying a panel that has been removed

- Rinse the panel in warm water; you may also use a mild detergent to clean it.
- Let the panel air dry.
- Reapply the panel to the white section on the front of the handle.

Attaching the brush head

- Align the brush head so the bristles point in the same direction as the front of the handle (Fig. 3).
- Firmly press the brush head onto the metal shaft until it stops.

Note: There is a small gap between the brush head and the handle.

Charging your Sonicare

- Plug the charger into a live electrical outlet.
- Place the handle on the charger (Fig. 4).
 - The battery charge indicator in the power on/off button flashes to indicate that the toothbrush is charging.
 - When the Sonicare is fully charged, the battery charge indicator stops flashing and remains on.

Note: If the battery charge of your Sonicare is low, you hear 3 beeps and the battery charge indicator flashes rapidly for 30 seconds after the brushing cycle.

Note: To keep the battery fully charged at all times, you may keep your Sonicare on the charger when not in use. When the battery is fully charged, the charger draws minimal energy from the electrical outlet. It takes at least 24 hours to fully charge the battery.

Tip: For extra convenience during travel, you can remove the travel charger and use it without the charger cover.

Using the Sonicare For Kids

Brushing instructions

For children under age 8, it is recommended that a parent or adult caregiver perform tooth brushing.

- Wet the bristles and apply a small amount of toothpaste.
- Place the toothbrush bristles against the teeth at a slight angle towards the gumline (Fig. 5).
- Press the power on/off button to turn on the Sonicare.
- Apply light pressure to maximize Sonicare's effectiveness and let the Sonicare toothbrush do the brushing for you.
- Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between the teeth. Continue this motion throughout the brushing cycle.

Note: To make sure you brush evenly throughout the mouth, divide the mouth into 4 sections using the KidPacer feature (see chapter 'Features').

- Begin brushing in section 1 (outside top teeth) a series of musical notes indicate it is time to move to section 2 (inside top teeth). At the next series of notes move to section 3 (outside bottom teeth). The next series of notes indicate it is time to move to section 4 (inside bottom teeth). The toothbrush will automatically turn off at the end of the brushing cycle (Fig. 6).

- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You may also brush your tongue, with the toothbrush turned on or off, as you prefer (Fig. 7).

Your Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

Brushing modes

The toothbrush automatically starts in the default high mode.

- Press the brushing mode button to toggle between modes (Fig. 8).

- The green mode light indicates the selected mode.
 - High mode: Full-strength mode designed for older/independent brushers.
 - Low mode: Reduced power mode used for training and introducing younger children to the sonic experience.

Features

KidTimer

- The KidTimer plays a congratulatory sound sequence at the end of the pre-determined brushing time. Additionally, the KidTimer automatically turns off the handle when the brushing cycle is complete.
- When in high mode, the KidTimer is set at 2 minutes.
- When in low mode the KidTimer is set at 1 minute and slowly increases over time. The progression of time helps young children reach the dentist recommended 2 minutes of brushing.

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush turns off.

KidPacer

The KidPacer uses a short series of musical notes to alert the brusher to move to the next brushing quadrant of the mouth. This ensures that the child gets a thorough, overall clean (see 'Brushing instructions').

Easy-start

- The Sonicare For Kids comes with the Easy-start feature activated.
- The Easy-start feature gently increases the power over the first 14 brushings.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.

Deactivating or activating the Easy-start feature

- Attach the brush head to the handle.

- Place the handle in the plugged-in charger.

To deactivate Easy-start:
Press and hold the power on/off button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

To activate Easy-start:
Press and hold the power on/off button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Sonicare's effectiveness in removing plaque.

Cleaning

Do not clean brush heads, handle, charger, charger cover or interchangeable panels in the dishwasher.

Toothbrush handle

- Remove the brush head and rinse the metal shaft area with warm water (Fig. 9).

Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- Use a damp cloth to wipe the entire surface of the handle.

Brush head

- Rinse the brush head and bristles after each use (Fig. 10).
- Remove the brush head from the handle and rinse the brush head connection at least once a week with warm water (Fig. 10).

Charger

- Unplug the charger before cleaning.
- Use a damp cloth to wipe the surface of the charger.

Storage

- If you are not going to use the Sonicare for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.
- You can store a brush head on the brushing head holder at the back of the charger.

Replacement

Brush head

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Use only Sonicare For Kids replacement brush heads.

Disposal

- The rechargeable battery inside your Sonicare cannot be replaced, but it is easy to remove for recycling. At the end of the toothbrush's life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below.
- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 11).
- The built-in rechargeable battery contains substances that may pollute the environment. Always remove the battery before you discard the appliance, and recycle or dispose of the battery according to local waste management requirements.
- If you have any recycling questions, contact your local waste management office.

Removing the rechargeable battery

Please note that this process is not reversible.

- To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Sonicare.

- Insert a screwdriver into the slot located in the bottom of the handle and turn counterclockwise until you see a gap between the bottom cap and the handle (Fig. 12).

- Insert the screwdriver into the gap and pry the bottom cap out of the handle (Fig. 13).

- Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 14).

- Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 15).

The rechargeable battery can now be recycled or disposed of and the rest of the product discarded appropriately.

Guarantee and service

Philips guarantees its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips expense provided that convincing proof of purchase in the qualifying period is provided.
Contact our Consumer Care Center at 1-800-682-7664 (North America), outside North America contact your local Philips Consumer Care Center.

WARRANTY EXCLUSIONS

What is not covered under warranty:





- Brush heads
- Interchangeable panels
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

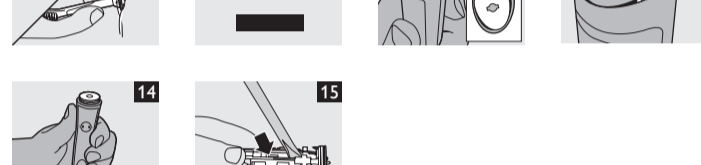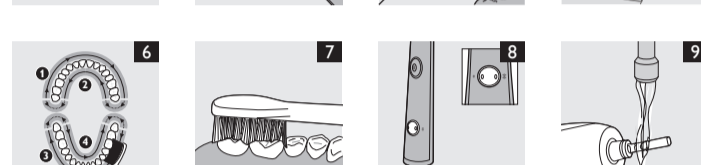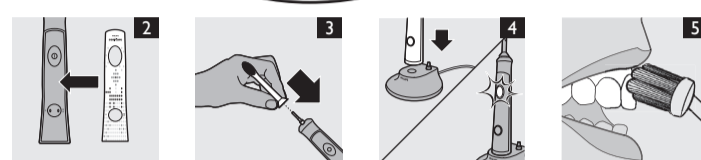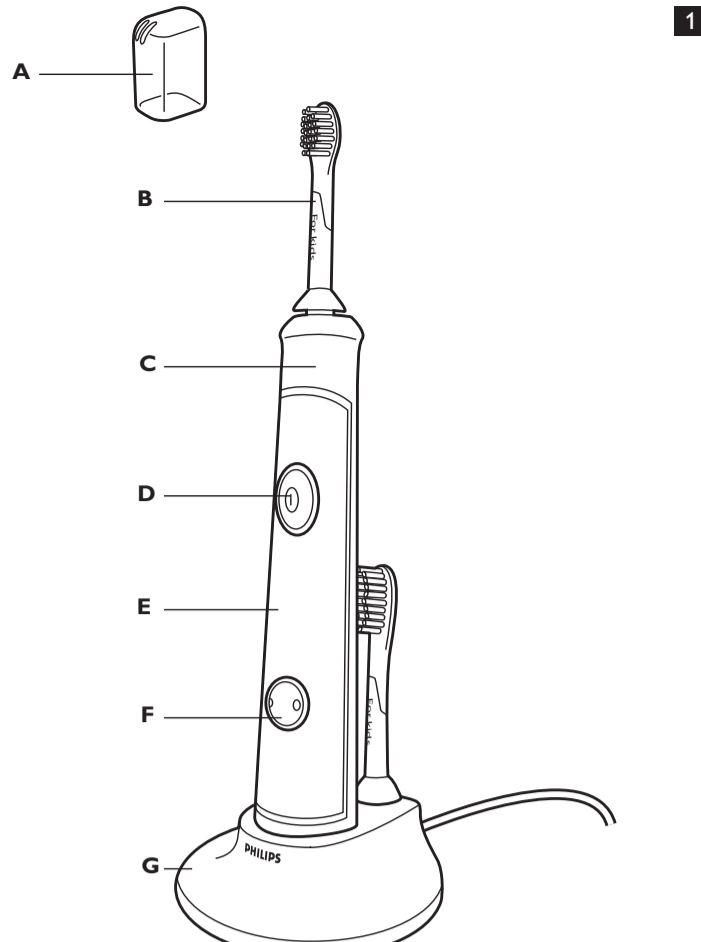IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

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©2009 Koninklijke Philips Electronics NV (KPEENV). All rights reserved. PHILIPS and the Philips shield are registered trademarks of KPEENV. Quadpacer and Sonicare are trademarks of Philips Oral Healthcare, Inc. and/or KPEENV.