# TANITA

Monitoring Your Health

# Digital Scale

Balance numérique  
Báscula (Pesa) Digital  
Balança digital

## Model HD-351

Operating Instructions / Mode d'emploi / Instrucciones para la Operación / Instruções operacionais

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### Federal Communications Commission and Canadian ICES Notice

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and Canadian ICES-003. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio or television technician for help.

### Modifications

The FCC requires the user to be notified that any changes or modifications made to this device that are not expressly approved by Tanita Corporation may void the user's authority to operate the equipment.

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### Features

- Previous weight recall.
- Digital accuracy and reliability.
- Classic simplicity fits any decor.
- Durable construction; non-skid surface.
- "Lo" message appears when batteries need replacing.


### Before Use

Select your preferred weight mode (lb, kg, or st. lb) using the Weight Mode Switch.

### Placing scale

The scale should be used on a hard even surface for accurate weighing. Not recommended for use on carpeted surfaces.

### Instructions for Use
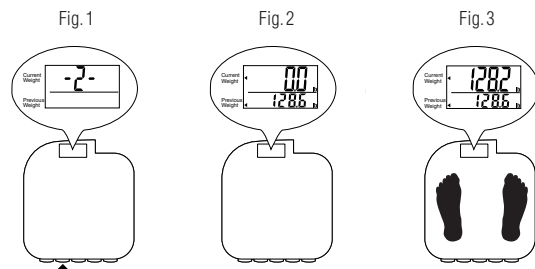
### To use scale

Tap the personal key with your toe.

The display will show personal number (fig.1)

Wait for the display to show "0.0" (fig.2), then step on the scale (fig.3).

Current reading in upper display and previous reading in lower display.

- Note :
- "OL" will appear on the display when the scale's weight capacity is exceeded.
  - "0.0" disappears if you do not step on the scale within 60 seconds.
  - The unit will shut off if you tap the Personal key 5 except while measuring.

### Precautions

To recalibrate the scale, gently tap the Personal Key with your toe (see Fig.4). Do not place your hand on the scale while recalibrating (see Fig.5); if you do, "-0.2" or "-0.4" will be displayed. Using your toe, tap the Personal Key again to display "0.0."

### Care and Maintenance

- Place the scale in an area free from heating equipment, high humidity, or extreme temperature changes.
- Never submerge the scale in water or spill chemical liquids on it. To clean, use a dampened cloth and mild detergent.
- Avoid excessive impact or vibration to the unit.

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