Oral-B®

Professional Care™

powered
by

BRAUN

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English 11

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Oral-B® Professional Care™ has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. However, when using electrical products, basic safety precautions should always be followed. Please read carefully before use.

Important

- Periodically check the cord for damage. A damaged or non-functioning appliance shall no longer be used. Take it to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age three.
- If the product is dropped, the brush-head should be replaced before the next use even if no damage is visible.

Description

- a Oral-B FlexiSoft brushhead
- b On/off button
- c Handle
- d Charge indicator
- e Charging unit
- f Brushhead compartment
- g Wall holder
- h Range of brushheads and accessories

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and is designed for use in the bathroom.

- Plug the charging unit (e) into an electrical outlet and place the toothbrush handle (c) on the charging unit.
- The **green** charge indicator (d) light flashes while the handle is being charged. If you then use your brush and place it back on the charging unit, the green charging light will again flash to let you know it is charging.


- When the handle is fully charged, the green charge indicator light shines continuously. A full charge takes at least 16 hours and enables approximately two weeks of regular brushing (twice a day, 2 minutes).
- The red charge indicator (d) light starts blinking when the rechargeable battery is low. Now it is time to recharge your toothbrush.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle by regular use at least every 6 months.

Specifications

Refer to the bottom of the charging unit for voltage specifications.

Using your toothbrush

Attach one of the brushheads (1–3) to the handle. They can be used with any kind of toothpaste. To avoid splashing, guide the brushhead to your teeth before switching the appliance on.

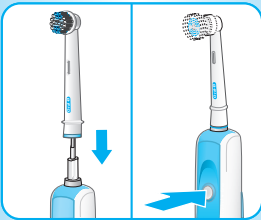D 19.545.2

«Clean»
«Sensitive»

D 19.585.3 X

«Clean»
«Sensitive»
«Polish»

Brushing modes

Your toothbrush offers different brushing modes to best meet your individual needs:

- «Clean» – Exceptional mouth cleaning.
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas.
- «Polish» – Exceptional polishing and natural whitening for occasional or everyday use.

How to switch between modes:

- Push the on/off button once, and your toothbrush automatically starts in the «Clean» mode.
- To shift from the «Clean» mode to the «Sensitive» mode, push the on/off button once more.
- If your toothbrush has a «Polish» mode, pushing the on/off button again will bring you from «Sensitive» to «Polish» mode.
- To turn your brush off from the «Polish» mode, push the on/off button once more.
- If you wish to turn your toothbrush off in the «Clean» or «Sensitive» mode, push and hold down the on/off button until the motor stops.

Outstanding Plaque Removal


For clean teeth and strong and healthy gums, **Oral-B Professional Care** gently cleans in areas where plaque bacteria can collect. It provides for an outstanding plaque removal even in hard-to-reach areas and is gentle on teeth and gums.


Brushing technique:

- Guide the brushhead slowly from tooth to tooth.
- Hold the brushhead in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

Pressure Sensor

For optimal brushing results **Oral-B Professional Care** has a pressure control feature to help you brush gently. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brushhead will

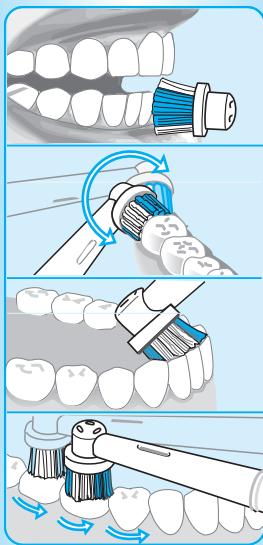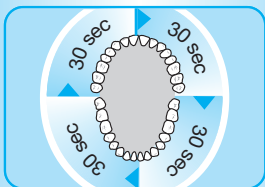continue but the pulsation will stop. In addition to feeling it, you can also hear a different sound.

Professional Timer

A short stuttering sound at 30-second intervals, which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a thorough consistent clean.

2-Minute Timer

A long stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.

Both timers memorize the elapsed brushing time, even when the handle is briefly switched off during brushing. Only when pausing longer than 30 sec. the timer resets to the initial starting point.


In the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Sensitive» mode.

Individual brushheads and accessories

The **Oral-B Professional Care** offers you a variety of brushheads and accessories, allowing you to select those which best meet your individual needs. They are all designed for daily use:

1 Oral-B FlexiSoft brushhead

The Oral-B FlexiSoft brushhead features outstanding oscillating-rotating action to effectively clean your teeth, clean along the gumline and in hard-to-reach areas such as behind teeth. The compact round brushhead surrounds each tooth for a thorough

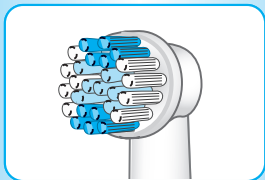
clean. The green FlexiSoft bristles bend, enabling both gentle brushing and deep reach of the blue Interdental Tips between teeth. Apply the brushing technique as described in the section «Outstanding Plaque Removal».

2 Oral-B Extra Soft brushhead

(not with all models)


The Oral-B Extra Soft brushhead provides thorough cleaning for sensitive teeth and gums. It has an extra soft bristle texture that is gentle on teeth and gums. Apply the same brushing technique as described in the section «Outstanding Plaque Removal».

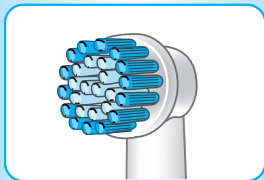
3 Oral-B Pro Bright brushhead

(not with all models)

The Oral-B Pro Bright brushhead is specifically designed to polish and naturally whiten your teeth by removing stains that can come from coffee, tea and tobacco. It features the oscillating-rotating head action of the Oral-B FlexiSoft brushhead to effectively clean teeth, combined with the special polishing cup for natural brightening. Apply the same brushing technique as described in the section «Outstanding Plaque Removal». We do not recommend using the Oral-B Pro Bright brushhead with braces. Children under 12 should not use the Oral-B Pro Bright brushhead.

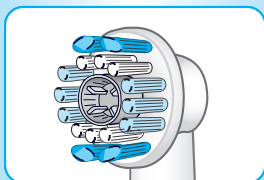
All three brushheads feature blue INDICATOR® bristles to help you monitor brushhead replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, signaling it is time to change the brushhead.

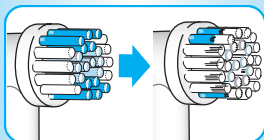
4 Interspace brushhead

(not with all models)

The Interspace brushhead is designed for cleaning between teeth as well as around bridges, crowns and implants. To use, position the brushhead at a

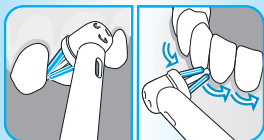45 degree angle to the gumline. Slowly guide the brushhead along the gumline, pausing between teeth for excellent results.

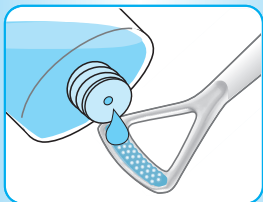5 Tongue Freshener


(not with all models)

The Tongue Freshener helps freshen breath. It is designed to help remove oral malodor and can be used everyday – before or after brushing. You can use mouthwash with the Tongue Freshener. Attach it to the handle, turn the appliance on and drag it across the tongue from back to front. We recommend using the «Sensitive» mode for enhanced control during use. Rinse after use. Replace the Tongue Freshener when the material becomes worn.

6 Interdental Cleaner

(not with all models)

The Interdental Cleaner effectively cleans between the teeth and helps keep your gums healthy. It should be used before brushing. Break one pick cleaner off the base (7) and insert it into the slot in the neck of the white pick holder. Attach the Interdental Cleaner to the handle and turn the appliance on. Depending on which part of your mouth you are cleaning, the pick holder should be adjusted by turning it, so that the blunt side always faces your gums. Guide the cleaner between your teeth and clean systematically from space to space. We recommend using the «Sensitive» mode for enhanced control during use. A pick cleaner can be used approximately two to three times until it becomes too soft to effectively fit and clean between the teeth. Replacement cleaners are available separately.

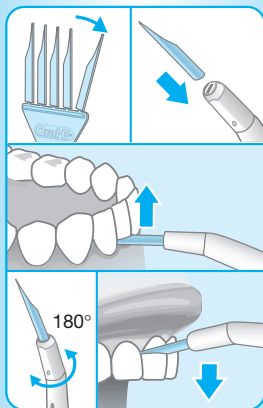Oral-B Professional Care replacement brushheads and accessories are available in stores that sell Oral-B products.

Cleaning recommendations


- After brushing, rinse your brushhead or accessory thoroughly under running water for several seconds with the handle switched on.
- Switch off the handle and remove the brushhead or accessory from the handle. Clean both parts separately under running water, then wipe them dry before you reassemble the brush and put the handle on the charging unit.
- From time to time, also clean the charging unit and the brushhead compartment with a damp cloth.
- For cleaning purposes, the brushhead compartment is removable.

Subject to change without notice.

This product conforms to the European Directives EMC 2004/95/EC and Low Voltage 2006/95/EC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

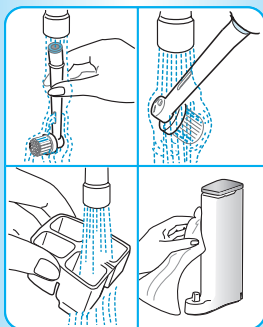Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brushheads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

Free Manuals Download Website

<http://myh66.com>

<http://usermanuals.us>

<http://www.somanuals.com>

<http://www.4manuals.cc>

<http://www.manual-lib.com>

<http://www.404manual.com>

<http://www.luxmanual.com>

<http://aubethermostatmanual.com>

Golf course search by state

<http://golfingnear.com>

Email search by domain

<http://emailbydomain.com>

Auto manuals search

<http://auto.somanuals.com>

TV manuals search

<http://tv.somanuals.com>