Oral-B®

ProfessionalCare™

8500 DLX OxyJet® Center

powered
by **BRAUN**

한국어	4
English	12

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The **Oral-B® ProfessionalCare™ 8500 DLX OxyJet® Center** is a dentist-inspired oral care system for cleaning, maintaining health teeth and gums, whitening and freshening.

The **Oral-B® ProfessionalCare™ 8500 DLX** toothbrush delivers outstanding performance:

- Unique 3D action for ultimate cleaning: 40.000 in-and-out pulsations gently loosen plaque, 8.800 side-to-side oscillations sweep plaque away
- Compact cup shaped brushhead: Surrounds each tooth to clean completely between and behind teeth
- Innovative accessories and brushheads to clean and polish teeth as well as freshen your tongue and mouth
- Customised Speed Control which offers a broad variety of speed levels to meet your individual brushing needs
- Premium handle design and ergonomics for greater comfort and control

The **Oral-B® ProfessionalCare™ OxyJet®** irrigator mixes air into water in the form of micro-bubbles. These micro-bubbles in the water stream attack plaque bacteria. At the same time, the irrigator removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

- The appliance must be mounted in a way that it cannot fall into water. The plug should always be accessible.
- Periodically check the cord for damage. A damaged or non-functioning appliance shall no longer be used. Take it to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age three.
- If the product is dropped, the brushhead should be replaced before the next use even if no damage is visible.

Description

Toothbrush

- ① Charger base
- ② Charge indicator
- ③ Toothbrush handle (rechargeable)
- ④ On/off switch
- ⑤ Brushhead (see «Brushheads and accessories»)
- ⑥ Brushhead compartment

Irrigator

- ⑦ Jet compartment
- ⑧ Irrigator jet
- ⑨ Jet switch
- ⑩ Jet release button
- ⑪ On/off switch
- ⑫ Water pressure control dial
- ⑬ Handle switch (water stop and go)
- ⑭ Irrigator handle
- ⑮ Water container (600 ml)
- ⑯ Cord storage

Specifications

Power supply: see bottom of housing

Power intake: 24 W

Using your Oral-B® ProfessionalCare™ 8500 DLX toothbrush

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and can be used in the bathroom without hesitation.

- Plug the appliance into an electrical outlet and place the toothbrush handle ③ on the charger base ①.
- The blue charge indicator ② flashes while the handle is being charged. If you then use your brush and place it back on the charger, the charging light will again flash to let you know it is charging.
- When the handle is fully charged, the charge indicator shines continuously. A full charge takes at least 16 hours and enables up to two weeks of regular brushing.
- When the rechargeable battery is running low, your toothbrush will slow down, then stop.
- For everyday use, the toothbrush handle can be stored on the charger base to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, refrain from placing the handle on the charger base and fully discharge it through regular use at least every 6 months.



Customised Speed Control

Oral-B® ProfessionalCare™ 8500 DLX offers a variable range of speeds to best meet your individual needs. High speed is ideal for tooth and gumline cleaning. For cleaning sensitive teeth and gum tissue as well as precision use of the accessories, we recommend using a lower speed.


- Press the on/off button once to operate at maximum speed.
- Continue to hold the on/off button to move from the highest speed to the Customised Speed Control mode.
- As long as you keep the on/off button pressed, the handle will continuously run through all speed levels.
- Release the on/off button to brush at your preferred speed level.
- Push the on/off button once more to turn the toothbrush off.
- If you then turn it on again you will automatically come back to the maximum speed level.



outsides

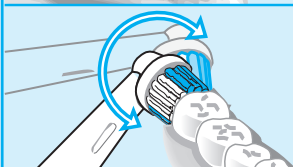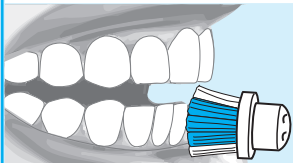insides

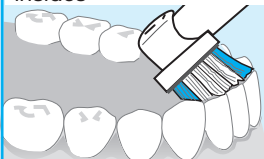tops

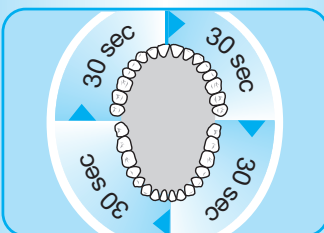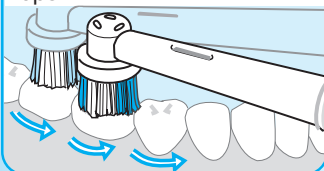Brushing technique

- Attach one of the brushheads ⑤ to the handle. They can be used with any kind of toothpaste.
- To avoid splashing, guide the brushhead to your teeth before switching the appliance on.
- Guide the brushhead slowly from tooth to tooth.
- Hold the brushhead in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

Pressure Sensor

For optimal brushing results **Oral-B® ProfessionalCare™ 8500 DLX** has a pressure control feature to ensure that the right degree of pressure is applied. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brushhead will continue but the pulsation will stop. In addition to feeling it, you can also hear the difference.

Professional Timer

A short stuttering sound at 30-second intervals, which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a better overall clean.

2-Minute Timer

A long stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.


During the first few days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist.

Brushheads and accessories

The **Oral-B® ProfessionalCare™ 8500 DLX** offers you a variety of brushheads and accessories to meet your individual needs. They are all designed for daily use:

I FlexiSoft® brushhead


The FlexiSoft® brushhead features our outstanding oscillating-pulsating action to effectively clean your teeth, clean along the gumline and in hard-to-reach areas such as behind teeth. The green FlexiSoft® bristles bend, enabling both gentle brushing and deeper reach of the blue Interdental Tips between teeth.

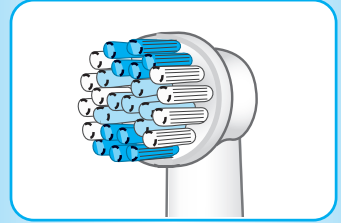
II ProBright™ brushhead

The ProBright™ brushhead is specifically designed to polish and naturally whiten your teeth by removing coffee, tea and tobacco stains. We do not recommend using the ProBright™ brushhead with braces.

Apply the same brushing technique as described in the section «Brushing technique» for both brushheads.

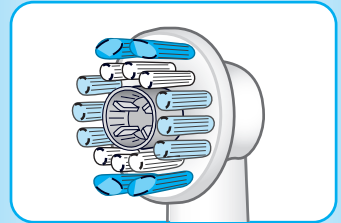
Both brushheads feature blue INDICATOR® bristles to help you monitor brushhead replacement need. With thorough brushing, twice a day for two minutes, the blue colour will fade halfway within 3 months, signaling it is time to change the brushhead.

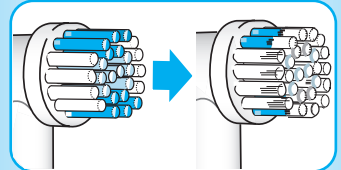III Tongue Freshener

The Tongue Freshener freshens breath by removing odour-causing bacteria from the tongue. It can be used everyday – before or after brushing. You can also use mouthwash with the Tongue Freshener.

Attach it to the handle, turn the appliance on and drag it across the tongue from back to front. We recommend a low speed setting for enhanced control during use. Rinse after use. Replace the Tongue Freshener when materials become worn.

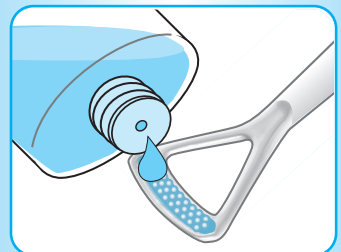IV Interdental Cleaner



The Interdental Cleaner effectively cleans between the teeth and helps improve your gum health. It should be used before brushing.

Break one pick cleaner off the base and insert it into the slot in the neck of the white pick holder. Then attach the Interdental Cleaner to the handle. Depending on which part of your mouth you are cleaning, the pick holder should be adjusted by turning it, so that the blunt side always faces your gums. Turn the appliance on and guide the cleaner between your teeth and clean systematically from space to space. We recommend a low speed setting for enhanced control during use. A pick cleaner can be used approximately two to three times until it becomes too soft to effectively fit and clean between the teeth. Replacement cleaners are available separately.

Oral-B® ProfessionalCare™ 8500 DLX replacement brushheads and accessories are available in stores that sell Oral-B products.

After use

- Rinse your brushhead or accessory thoroughly under running water for several seconds with the handle switched on.
- Switch off the handle and remove the brushhead or accessory from the handle. Clean all parts separately under running water, then wipe them dry before you reassemble the brush and put the handle on the charger base.

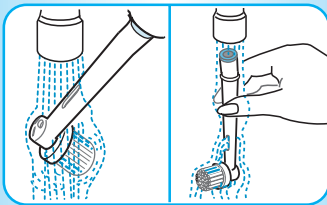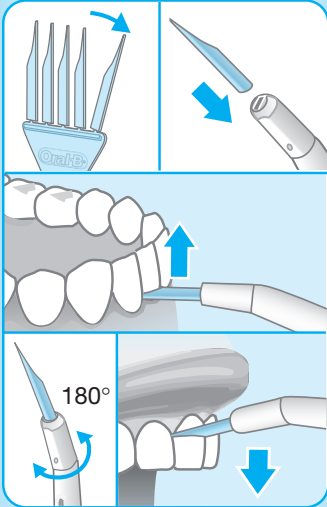Using your Oral-B® ProfessionalCare™ OxyJet® irrigator

Before using for the first time, fill the water container ⑮ and allow the water to run through completely by holding the irrigator handle ⑯ over a washbasin.

1. Attach one of the irrigator jets ⑰ to the handle until it locks into place. Fill the water container with water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch ⑱. (When doing so, the handle switch should be in the lower position = water stop).
3. Select the water pressure with the pressure control dial ⑲ (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch ⑳ up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.







You can choose between two different jet functions by sliding the jet switch ㉑ down or up (before doing so, turn off the water stream by sliding the handle switch ⑳ down):


Rotating stream: for general cleaning and gum massage (jet switch in lower position). The rotation of the water stream and its micro-bubbles is designed for easy penetration.

Mono stream: straight water stream for site-specific cleaning (jet switch in upper position).

Max. operating time: 15 minutes
Cooling off time: 2 hours

To ensure the jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

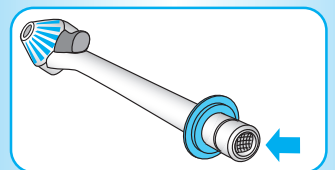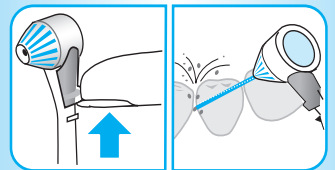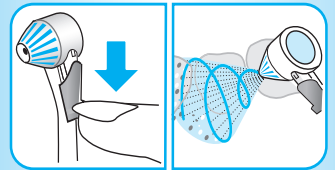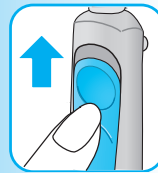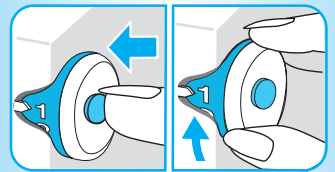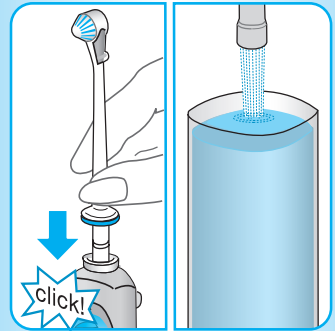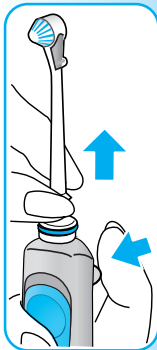After use

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended the use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator. It will automatically turn off after 10 minutes if you forget to turn it off manually. In this case, to turn it on again, you have to press the on/off switch twice.

To take off the irrigator jet, press the jet release button [®]. The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brushheads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Customer Service Centre.

Free Manuals Download Website

<http://myh66.com>

<http://usermanuals.us>

<http://www.somanuals.com>

<http://www.4manuals.cc>

<http://www.manual-lib.com>

<http://www.404manual.com>

<http://www.luxmanual.com>

<http://aubethermostatmanual.com>

Golf course search by state

<http://golfingnear.com>

Email search by domain

<http://emailbydomain.com>

Auto manuals search

<http://auto.somanuals.com>

TV manuals search

<http://tv.somanuals.com>