BRAUN
Oral-B

Vitality

D 12013
D 12013W

Type 3709

3-709-377/00/IV-07
Japan/GB

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English

Please read the use instructions carefully before first use.

Important

- Periodically check the cord for damage. A damaged or non-functioning unit shall no longer be used. Take the charging unit (E) to a Braun Repair Centre.
- This product is not intended for use by children under age three.
- If the product is dropped, the brushhead should be replaced before the next use even if no damage is visible.

Description

- A Brushhead (depending on model) – interchangeable
- B On switch
- C Off switch
- D Handle
- E Charging unit

Charging

- Plug the charging unit (E) into an electrical outlet. Place the handle on the charging unit. Handle must be switched off.
- A full charge takes 16 hours and will provide an operation time of approximately 20 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, you must discharge it at least once every 6 months. To do this, unplug the charging unit and allow the battery to run down with normal use.

Using the toothbrush

- Wet brushhead and apply toothpaste.
- Place brushhead on teeth, then switch on.
- Slowly guide the brushhead from tooth to tooth, first the outsides, then the insides, finally the chewing surfaces. Do not press too hard or scrub.

- After two minutes (brushing time recommended by dentists) the timer (D 12513 model only) will signal by a short stutter.
- During the first days of use, gums may bleed slightly. Should bleeding persist after 2 weeks, contact your dentist.
- Replace brushhead every 3 months, when bristles are worn (1c) or when colour on Indicator bristles (1b) has faded half-way down the bristles.

Cleaning


After use, rinse the brushhead under running water (2).

Remove the brushhead and rinse brushhead and handle separately (3), then dry.

Environmental notice

This product contains nickel-cadmium rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Repair Centre.

However, if you prefer to dispose of the battery yourself, open the handle as shown below:

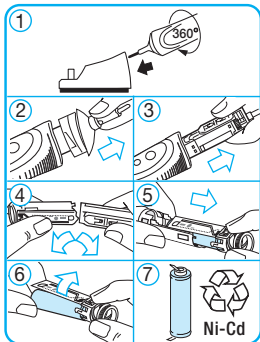Caution! Opening the handle will destroy the appliance and invalidate the guarantee.

Subject to change without notice.

Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brushheads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Repair Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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