# Professional Care

D 16.565

D 16.525

D 16.525 H

D 16.523

D 16.513







powered by **BRAUN**

Deutsch	4	(DE) (AT)	00 800 27 28 64 63
English	7	(CH)	08 44 - 88 40 10
Français	12	(GB)	0 800 783 7010
Español	16	(IE)	1 800 509 448
Português	19	(FR)	0 800 944 803 (appel gratuit depuis un poste fixe)
Italiano	22		
Nederlands	25	(BE)	0 800 14 592
Dansk / Norsk	28	(LU)	4 00 50 51
Svenska	31	(ES)	901 11 61 84
Suomi	34	(PT)	808 20 00 33
Ελληνικά	37	(IT)	800 124 600
		(NL)	0800-4 45 53 88
		(DK)	70 15 00 13
		(NO)	22 72 88 10
		(SE)	020 - 21 33 21
		(FI)	020 377 877
<b>Internet:</b> www.oralb.com www.braun.com www.service.braun.com		(GR)	800 801 3457
		(AU)	1 800 641 820
		(NZ)	0 800 108 909
Braun GmbH Frankfurter Straße 145 61476 Kronberg/Germany		(HK)	852-25249377 (Audio Supplies Company Ltd.)

## Type 4729
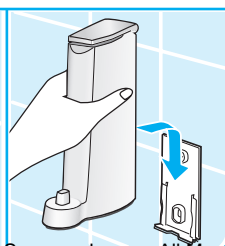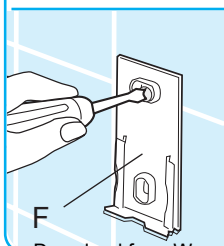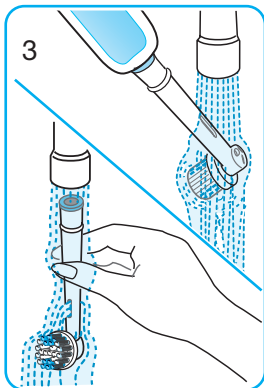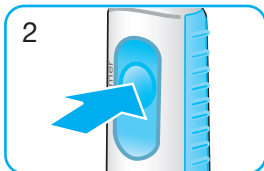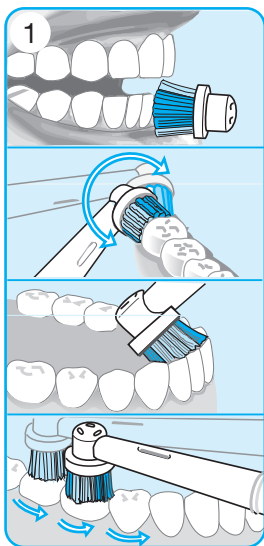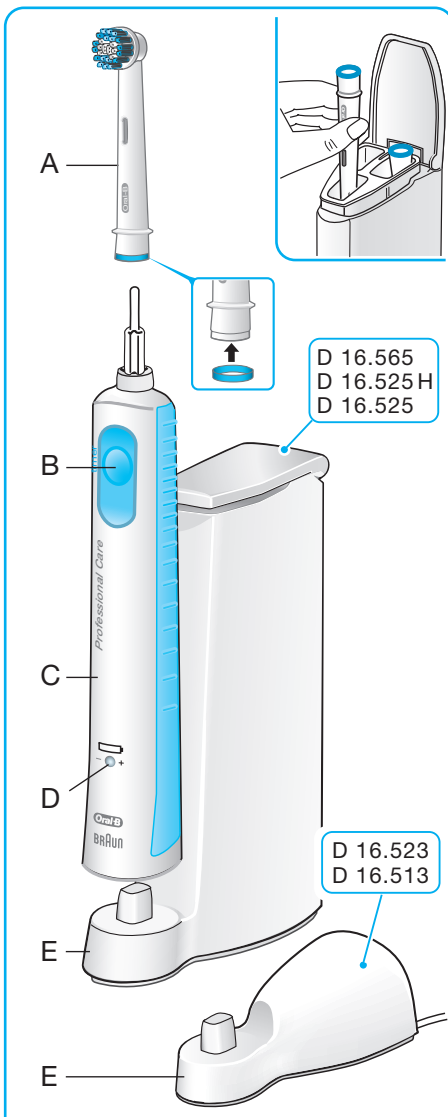
## Deutsch

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CH: 0844 - 88 40 10

## English

Please read the usage instructions carefully before using first time.

### Important

- Periodically check the cord for damage. If the cord is damaged, take the charging unit to an Oral-B Braun Service Centre. A damaged or non-functioning unit must no longer be used.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.
- Children should be supervised to ensure that they do not play with the appliance.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service Centre. Do not disassemble the product except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not insert any object into any opening of the appliance.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use accessories not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

### Description

- A Brush head
- B On/off button
- C Handle
- D Charge indicator
- E Charging unit

### Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (E) into an electrical outlet and place the toothbrush handle (C) on the charging unit.
- The green charge indicator light (D) flashes while the handle is being charged, once it is fully charged the green charge indicator light shines continuously for 5 seconds and then it turns off. A full charge takes at least 18 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).

- When the rechargeable battery is low, the red charge indicator (D) light starts flashing after turning off.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible. However, for environmental reasons, Oral-B recommends to unplug the charger until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

## Specifications

Power supply: see bottom of charging unit

Handle voltage: 1.2 V

## Using the toothbrush

You can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on the appliance. Guide the brush head from tooth to tooth, spending a few seconds on each tooth surface. Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces. Do not press too hard or scrub (1).

## Pressure Sensor

To promote optimal brushing your toothbrush has a pressure control feature installed. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition you will also hear a different sound while brushing.

## Professional Timer

A short stuttering sound at 30-second intervals reminds you to brush equally all four quadrants of your mouth. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. The elapsed brushing time is memorised when the handle is briefly switched off during brushing. When pausing for longer than 30 seconds the timer resets.

In the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or dental hygienist.

## Cleaning

After use, rinse the brush head for several seconds under running water with the handle switched on. Then switch off the handle and remove the brush head. Rinse both parts separately under running water (3) and wipe them dry. From time to time, also clean the charging unit with a damp cloth (4).

Subject to change without notice.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

## For UK/IRL only

### Warranty

We grant 2 years warranty on the product commencing on the date of purchase. Within the warranty period, we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This warranty extends to every country where this appliance is supplied by Braun or its appointed distributor.

This warranty does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The warranty becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the warranty period, hand in or send the complete appliance with your original till receipt to an authorised Oral-B Braun Customer Service Centre.

#### UK:

FREEPOST RSHJ-KRSA-UHGE

Braun Product Service

Olympic Shaver Centre Ltd, PO BOX 9799, Calverton

Nottingham, NG14 6WN

#### ROI:

Braun Product Service

Olympic Shaver Centre Ltd, Dept AA1758, PO BOX 6151

FREEPOST, Dublin 2

Or please call the Oral-B Customer Service Line if you have any queries:

### Oral-B Customer Services

UK (free phone): 0 800 783 7010

IRL (free phone): 1 800 509 448

This warranty in no way affects your statutory rights!

## For AUS/NZ only

### Warranty

We grant a 24 Month Warranty on this appliance commencing on the date of purchase.

**This warranty applies in all States and Territories of Australia and in New Zealand and is in addition to any mandatory statutory obligations imposed on Procter & Gamble Australia Pty Ltd or Procter & Gamble Distributing New Zealand, its distributors and its manufacturers. This express warranty does not purport to exclude, restrict or modify any such mandatory statutory obligations.**

This warranty only covers defects in the appliance resulting from faults in workmanship or materials. Therefore, if your appliance becomes faulty as a result of faults in workmanship or materials, it will be either exchanged by

your retailer, or repaired/exchanged by our Service Agent, on presentation of proof of purchase.

This warranty only applies to domestic or household use of this appliance and the warranty will only apply if the correct operating instructions as set out in this manual have been followed. The warranty on any appliance replaced under this warranty ends on the expiry of the warranty period that applied to the original purchase of appliance.

This warranty does not cover:

- A. Damage arising from improper use or operation on incorrect voltages.
- B. Breakages to replacement parts (e.g. brush heads, shaver foils, glass jars, etc.).
- C. Normal wear due to moving parts.
- D. Repairs undertaken by unauthorised service personnel or use of non genuine parts.
- E. Appliances that are outside the warranty period or are not faulty.
- F. Any failure to provide facilities for repair or spare parts more than 2 years after the date of purchase of appliance.

This appliance is not intended for use by young children or infirm persons without supervision. Young children should be supervised to ensure that they do not play with the appliance.

If the appliance is powered by a cord extension set or electrical portable outlet device these should be positioned so they are not subject to splashing or ingress of moisture.

To the extent permitted by the laws of the commonwealth of Australia and New Zealand and applicable laws of each State or Territory of Australia, and subject to the second paragraph of this warranty above, this warranty is given in exchange for any other rights you may have against Procter & Gamble Australia Pty Ltd or Procter & Gamble Distributing New Zealand or its distributors or manufacturers, whether at law, in tort (including negligence), in equity or under statute.

Please retain proof of purchase of this appliance. If you have any questions in relation to this warranty, please call our Consumer Service line (see below for numbers).

### **Repairs and service**

If you are sending your appliance for repair, please pack the appliance adequately and send it to your nearest Service Agent.

For service or replacement parts please refer to the authorised Service Agents listed overleaf.

Local Service Agents may change from time to time – in that event please contact the Consumer Service freecall number below for updated local Service Agent details.

### **Australian Consumer Service**

**(free call): 1 800 641 820**

### **New Zealand Consumer Service**

**(free call): 0 800 108 909**

N2820

## **Australian & New Zealand Service Agents**

### **VICTORIA & TASMANIA**

Statewide Services Pty Ltd  
(trading as J A Appliances)  
17-19 Hossack Avenue  
Coburg North, VIC 3058  
Tel. (03) 8371 9100  
Fax: (03) 9355 8644  
Email: jaappliances@statewideservices.com.au

### **QUEENSLAND**

Walters Import trading as  
Brisbane Appliance Service Centre  
449 Beaudesert Road  
Moorooka  
QLD 4105  
Tel. (07) 3277 0360  
Fax. (07) 3274 2909  
Email: basc@waltersimport.com.au

### **WESTERN AUSTRALIA**

Mark 1 Appliances  
Unit 1, 46 Abernethy Road  
Belmont, WA 6104  
Tel. (08) 9475 0911  
Fax: (08) 9475 0922  
Email: mark1app@iinet.net.au

### **NEW SOUTH WALES & ACT**

Faulder Appliance Services Pty Ltd  
Shop 4, 9 Evans Road,  
Rooty Hill, NSW 2766  
Tel. (02) 9625 3064  
Fax: (02) 9625 3037  
Email: fas16@bigpond.com

### **SOUTH AUSTRALIA & NT**

Adelaide Service Centre  
Shop B2,  
Parkholme Shopping Centre  
319 Oaklands Road, Parkholme  
Adelaide, SA 5043  
Tel. (08) 8357 5800  
Fax. (08) 8357 5833  
Email: ascckh@adam.com.au

### **NEW ZEALAND**

Key Service Ltd  
7D Echelon Place  
East Tamaki  
Manukau City  
Auckland  
New Zealand  
(PO Box 38569, Howick,  
Auckland, New Zealand)  
Tel. (09) 916 0960  
Fax: (09) 916 0970  
Email: key-service-ltd@xtra.co.nz

## **30-Days Money Back Guarantee**

Try Oral-B Professional Care starting from the day of purchase. If you are not satisfied with this product, receive a full refund. For a full refund, please send us the product along with your name, address, contact phone number and the original cash register receipt indicating purchase price and date to the address below.

UK: Oral B Money Back Guarantee  
Dept. PX415, PO Box 10, Diss, Norfolk, IP98 3SP

ROI: Oral B Money Back Guarantee  
Dept. PX415, PO Box 1760, Tallaght, Dublin 24

This money back guarantee in no way affect your statutory rights.

Australia: «Oral-B 30 Day Trial Offer»  
PO Box 150, Brookvale NSW 2100, Australia

New Zealand: «Oral-B 30 Day Trial Offer»  
PO Box 33083, Takapuna, Auckland, New Zealand

Your Oral-B Triumph must be sent in its original packaging and be postmarked no later than 30 days after the date of purchase. Please allow 6-8 weeks for handling. The shipment cost to return the product is the responsibility of the consumer.

Sorry: No refunds will be made for any Oral-B Professional Care product damaged by accident, neglect, unreasonable use or lost in the post.

Terms and conditions apply. For further details please contact our customer service.

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