
Oral-B®

**AdvancePower
WaterJet**

MD 31

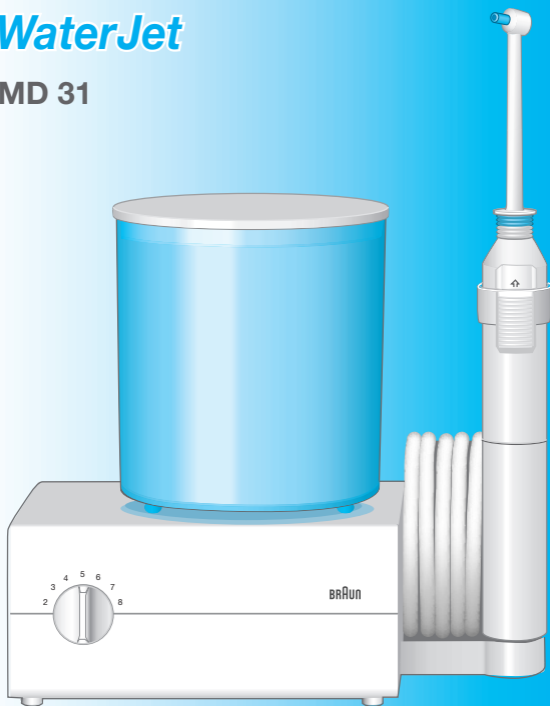powered
by

BRAUN

Type 4803










Deutsch	4, 22, 27
English	7, 22, 27
Français	10, 23, 27
Español	13, 24, 27
Italiano	16, 24, 27
Nederlands	19, 25, 27

Internet:

www.oralb.com
www.braun.com

- (D)** Haben Sie Fragen zu diesem Produkt?
Rufen Sie an (gebührenfrei):
00 800 27 28 64 63
- (A)**
- (CH)** **08 44 - 88 40 10**
- (GB)** If you should require further assistance, please call our Consumer Service line.
0 800 73 11 792
(free call)
- (IRL)** **Helpline**
1 800 509 448 (free call)
- (F)** Appelez le Service Consommateurs Oral-B :
0 810 365 855
(numéro azur)
- (B)** **02-711 92 11**
- (E)** Servicio al consumidor para España:
901 11 61 84
- (I)** Servizio consumatori:
02-667 86 23
- (NL)** Heeft u vragen over dit produkt?
Bel Braun/Oral-B Consumenten-infolijn:
0800-4 45 53 88

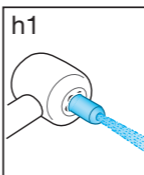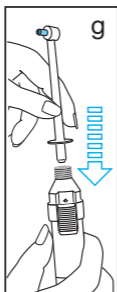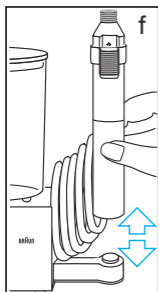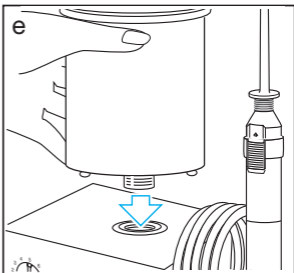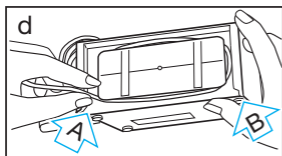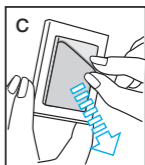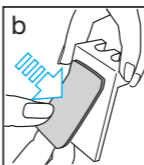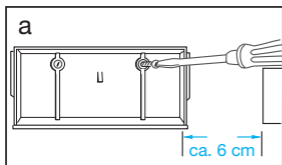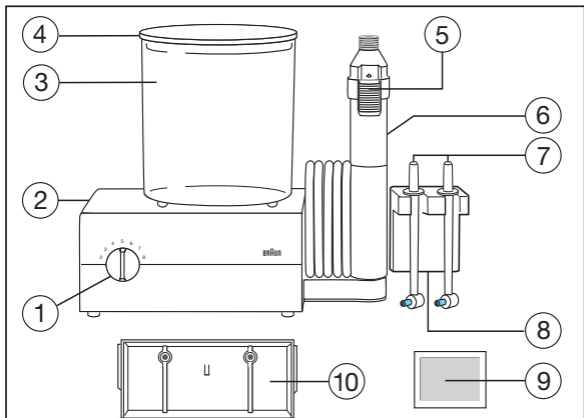Deutsch

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English

The Oral-B AdvancePower WaterJet removes remaining food particles from between your teeth, while it also effectively massages your gums.

Please read the use instructions carefully before first use.

Important

The appliance must be mounted in a way that it cannot fall into water.

The plug should always be accessible. Periodically check the cord for damage. A damaged or non-functioning unit should no longer be used. Take it to a Braun Service Centre.

Specifications

Max. operating time: 20 minutes

Cooling off time: 2 hours

Water capacity: approx. 500 cm³

Description

- ① Water impact control
- ② Motor part
- ③ Water container
- ④ Lid
- ⑤ Magnetic ring switch (removable)
- ⑥ Hand piece
- ⑦ 2 attachable jets with different colour markings
- ⑧ Wallholder for the jets
- ⑨ Double-faced adhesive strip for the wallholder
- ⑩ Wallholder for the unit

Wallholder

The appliance comes with a wallholder ⑩ for the unit ②. Use 30 x 3.5 mm semi-round screws with wall plugs, if necessary, to mount the holder to the wall (a). Slide the unit onto the holder from the top.

The wallholder ⑧ for the jets ⑦ can be mounted with the double-faced adhesive strip ⑨. Clean the tiles thoroughly with surgical or methylated spirit first. Then pull off the foil from one side of the adhesive strip and press this side against the back of the wallholder (b). Then pull off the foil from the other side of the strip (c) and press the wallholder firmly against the wall. Wait 12 hours before hanging up the jets to ensure maximum adhesion.

Switch

The magnetic ring switch ⑤ on the hand piece turns the WaterJet on and off:

Off: Switch is below the arrow

On: Slide switch over the arrow

The water impact control ① regulates the impulse pressure:
gentle = 2, strong = 8

The water impact is continuously variable and is electronically controlled.

Setting in operation

The cord compartment offers practical storage for the cord (d). You only have to use as much cord as you need to reach the socket. Guide the cord through either the left- or right hand opening of the cord compartment (A/B).

After connecting the unit to a socket, remove the water container ③ and fill it with lukewarm water. Then replace the water container (with the lid ④ on) firmly onto the unit (e). Before using the first time, allow the water to run through completely. Take the hand piece ⑥ out of the magnetic holder (f), hold it into the washbasin and turn the WaterJet on.

Using

Fill the water container ③ with water again and replace it firmly onto the unit ② (e). Attach one of the jets ⑦ onto the hand piece ⑥ and let it audibly click in place (g). Select the water pressure you want with the water impact control ①; it is best to start with the gentlest pressure (setting 2), especially if you tend to have bleeding gums or if children will be using the WaterJet.

Bend over the washbasin and guide the jet to your teeth.

Then turn the unit on by activating the switch ⑤. The WaterJet impulses will start in a few seconds. Keep your mouth opened slightly to allow the water to run into the washbasin.

Mono- and multi-jet irrigation

The jet nozzle offers both mono- and multi-jet irrigation (h):

Mono-jet irrigation: Jet nozzle in position (h1). Cleans teeth and the spaces between them, flushing out even stubborn left-over food particles.

Multi-jet irrigation: Jet nozzle in position (h2). Effectively massages gums to keep them strong and healthy.

Start with mono-jet irrigation. Direct the jet nozzle to your teeth and especially the spaces in between. Then massage your gums with multi-jet irrigation.

You can easily switch over from one setting to the other. Simply turn the WaterJet off with the switch ⑤ (to avoid splashing) and press the jet nozzle in or out.

Then turn the WaterJet on again with the switch ⑤.

After use

Do not turn the unit off before the water container is entirely empty. Bacteria could otherwise develop in the water if the unit is not used over a longer period of time.

After use, replace the hand piece ⑥ on the magnet holder (f).

Cleaning

Clean the motor part ② with a damp cloth only. Wipe the hand piece ⑥ dry. Cleaning it under the switch ⑤ presents no problems, because the switch can be pulled off completely towards the top.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

At the end of the product's useful life, please dispose of it at appropriate collection points provided in your country.

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English

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends

to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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Hong Kong

Audio Supplies Company,
Room 506,
St. George Building,
2 ICE House Street,
Hong Kong,
☎ (852) 2524 9377

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Tehran Bouran Company
145 Ghaem Magham Farahani
Tehran
☎ +982 18 31 27 66

Ireland (Republic of)

Gavin's Electronics,
83-84, Lower Camden Street,
Dublin 2,
☎ 1800 509 448

Israel

S. Schestowitz Ltd.,
8 Shacham Str.,
Tel-Aviv, 49517,
☎ 1 800 335 959

Italia

Servizio Consumatori Braun
Gillette Group Italy S.p.A.,
Via G. B. Pirelli, 18,
20124 Milano,
☎ 02 / 6678623

Jordan

Interbrands
Wadi Al-Sir, Adwiat alhekma
St. Mekadabi bld.
Amman
☎ +692 582 75 67

Kenya

Radbone-Clark Kenya Ltd.,
P.O. Box 40833,
Mombasa Road,
Nairobi,
☎ 2 82 12 76

Korea

Gillette Korea Ltd.
144-27 Samsung-dong,
Kangnam-ku,
Seoul, Korea,
☎ 080-920-6000

Kuwait

Union Trading Company,
Braun Service Center,
P.O. Box 28 Safat,
Safat Code 13001, Kuwait,
☎ 04 83 32 74

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Philippines
Gillette Philippines Inc.,
Corporate Corner Commerce
Avenues
20/F Tower 1, IL Corporate Centre
1770 Muntinlupa city
☎ 02-771071 02-06/-16

Poland
Gillette Group Poland
Sp. z o.o.
ul. Domaniewska 41,
02-672 Warszawa,
☎ 22 548 88 88

Portugal
Grupo Gillette Portugal, Lda.,
Braun Service,
Rua Tomás da Fonseca,
Torre G-2°A,
1600-209 Lisboa,
☎ 808 2 000 33

Réunion
Dindar Confort,
Rte du Gymnase,
Boite Postale 278,
97940 St. Clotilde,
☎ 026 92 32 03

Romania
Gillette Romania srl.
Calea Floreasca nr. 133-137
et 1, sect 1,
714011 Bucuresti
☎ 01-2319656

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Singapore (Republic of)
Beste (S) Pte. Ltd.,
No. 6 Tagore Drive,
03-04 Tagore Industrial Building,
Singapore 787623,
☎ (65) 6552 2422

Slovakia
Techno Servis Bratislava
Trhová ul.
841 02 Bratislava IV
☎ 02 /64 46 36 43

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South Africa (Republic of)
Fixnet After Sales Service,
159 Queen Street, Kensington South,
P.O. Box 751770,
Johannesburg 2094,
Kensington South,
☎ 11 615 6765

St. Maarten
Rupchand Sons n.v. (ram's),
P.O. Box 79,
Philipsburg,
Netherlands Antilles,
☎ 52 29 31

St. Thomas
Boolchand's Ltd.,
31 Main Street,
P.O. Box 5667,
00803 St. Thomas,
US Virgin Islands,
☎ 340 776 0302

Suomi
Gillette Group Finland Oy,
P.O. Box 9,
Niittykatu 8, PL 9, 02200 Espoo,
☎ 09-45 28 71

Sverige

Gillette Group Sverige AB,
Dept. Sweden,
Stockholm Gillette
Räsundavägen 12,
Box 702,
16927 Solna,
☎ 020-21 33 21

Syria

Ahmed Hadaya Company
Hadaya building
Ain Keresh
Unisyria, P.O. Box 35002,
Damascus,
☎ 963 011-231433

Taiwan

Audio & Electr. Supplies Ltd.,
Brothers Bldg., 10th Floor,
85 Chung Shan N Rd., Sec. 1,
Taipei (104),
☎ (886) 02 2523 3283

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Yemen (Republic of)

Saba Stores for Trading,
26th September Street,
P.O. Box 5278,
Taiz,
☎ 4-25 23 88

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