

Note: The most current version of these instructions can be found online at www.grahamfield.com

- ⚠ **WARNING: Important! Read and understand these instructions before assembling the Versacart. If you do not understand any part of these warnings, cautions or instructions, contact a healthcare professional for direction in the use of this product.**
- ⚠ **WARNING: If components are damaged or missing, contact your dealer immediately. DO NOT use substitute parts. Use only Lumex® replacement parts. Non-Lumex® replacement parts could cause personal injury and damage to the Versacart.**
- ⚠ **WARNING: GF Health Products, Inc. assumes no responsibility for any damage or injury caused by improper assembly or use of this product.**


ASSEMBLY

Unpacking

1. Check for any obvious damage to the carton or its contents. If damage is evident, notify the carrier or your authorized Graham-Field Dealer.
2. Remove all loose packing from the carton.
3. Carefully remove all the components from the carton.

Versacart Assembly

1. Position Versacart frame in upright position on casters as shown at right. It may be necessary to first pull Versacart lower legs apart so it can stand upright.
- ⚠ **WARNING: Take care not to be pinched by moving parts while opening Versacart.**
2. Open handles and push down on upper and lower locking bars to secure Versacart frame in open, locked position as shown at right.

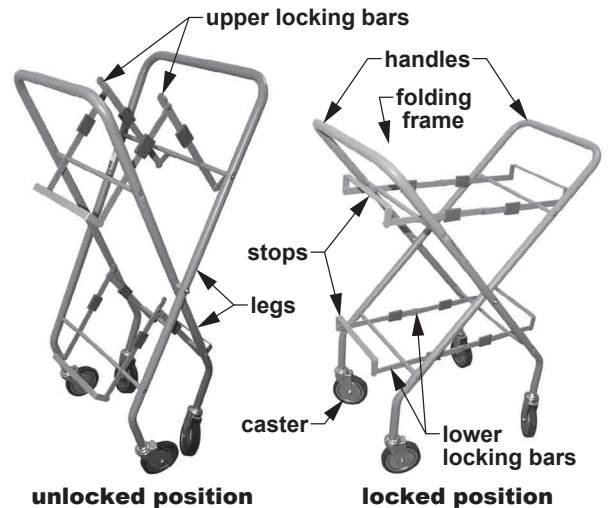

7010A Versacart assembly

BASKET INSTALLATION

Lower basket: Slide basket between legs onto lower locking bars. Position basket securely between stops as shown at lower right.

Upper basket: Position basket on upper locking bars securely between stops as shown at lower right.

- ⚠ **WARNING: Ensure that Versacart and Wire Baskets are assembled as described before use.**

OPERATION

- ⚠ **WARNING: The Versacart's maximum weight capacity is 200 lb. Do not exceed this limit.**



7010A Versacart assembled

MAINTENANCE

Check the fit and tightness frequently of all fasteners to ensure that the Versacart is securely assembled.

Cleaning

1. To clean: Gently wipe the Versacart with a soft, clean, damp cloth.

▲ **CAUTION: Do not use cleansers that could damage the Versacart's finish.**

2. To wash more thoroughly: Clean with mild soap and water.

3. Dry Versacart with a clean, soft cloth.

SPECIFICATIONS

Model	7010A Versacart
Support Width	24"
Support Depth	9.9"
Height	41"
Base Width	24"
Base Depth	16.5"
Casters	5" swivel
Weight	29 lb
Weight Capacity	200 lb

Model	7020A Wire Basket, Small
Size (W x L x H)	14" x 24" x 6"

Model	7021A Wire Basket, Large
Size (W x L x H)	14" x 24" x 12"

WARRANTY

GF Health Products, Inc. ("Graham-Field") warrants the Lumex® 7010A Versacart Folding Utility Cart and its components to be free from defects in workmanship and materials as listed below:

Versacart Frame, 7020A Basket, 7021A Basket: Three-Year Limited Warranty

Casters: One-Year Limited Warranty

If a product is deemed to be under warranty, GF Health Products, Inc. shall provide, at its option, (1) replacement of any defective part or product or (2) a credit of the original selling price made to GF Health Products, Inc.'s initial customer. The warranty does not include any labor charges incurred in replacement part(s) installation or any associated freight or shipping charges to GF Health Products, Inc. This warranty does not cover product failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and tear and need periodic replacement.

The warranties contained herein contain all the representations and warranties with respect to the subject matter of this document, and supersede all prior negotiations, agreements and understandings with respect thereto. The recipient of this document hereby acknowledges and represents that it has not relied on any representation, assertion, guarantee, warranty, collateral contract or other assurance, except those set out in this document.



GF Health Products, Inc. is not responsible for typographical errors. Packaging, warranties, products and specifications are subject to change without notice. Graham-Field and Lumex are registered trademarks of GF Health Products, Inc. © 2010 GF Health Products, Inc., 2935 Northeast Parkway, Atlanta, GA 30360, tel 800-347-5678, fax 800-721-0601

Free Manuals Download Website

<http://myh66.com>

<http://usermanuals.us>

<http://www.somanuals.com>

<http://www.4manuals.cc>

<http://www.manual-lib.com>

<http://www.404manual.com>

<http://www.luxmanual.com>

<http://aubethermostatmanual.com>

Golf course search by state

<http://golfingnear.com>

Email search by domain

<http://emailbydomain.com>

Auto manuals search

<http://auto.somanuals.com>

TV manuals search

<http://tv.somanuals.com>